



UltraFlex™ 7.5

Wire-Reinforced for Kink Resistance

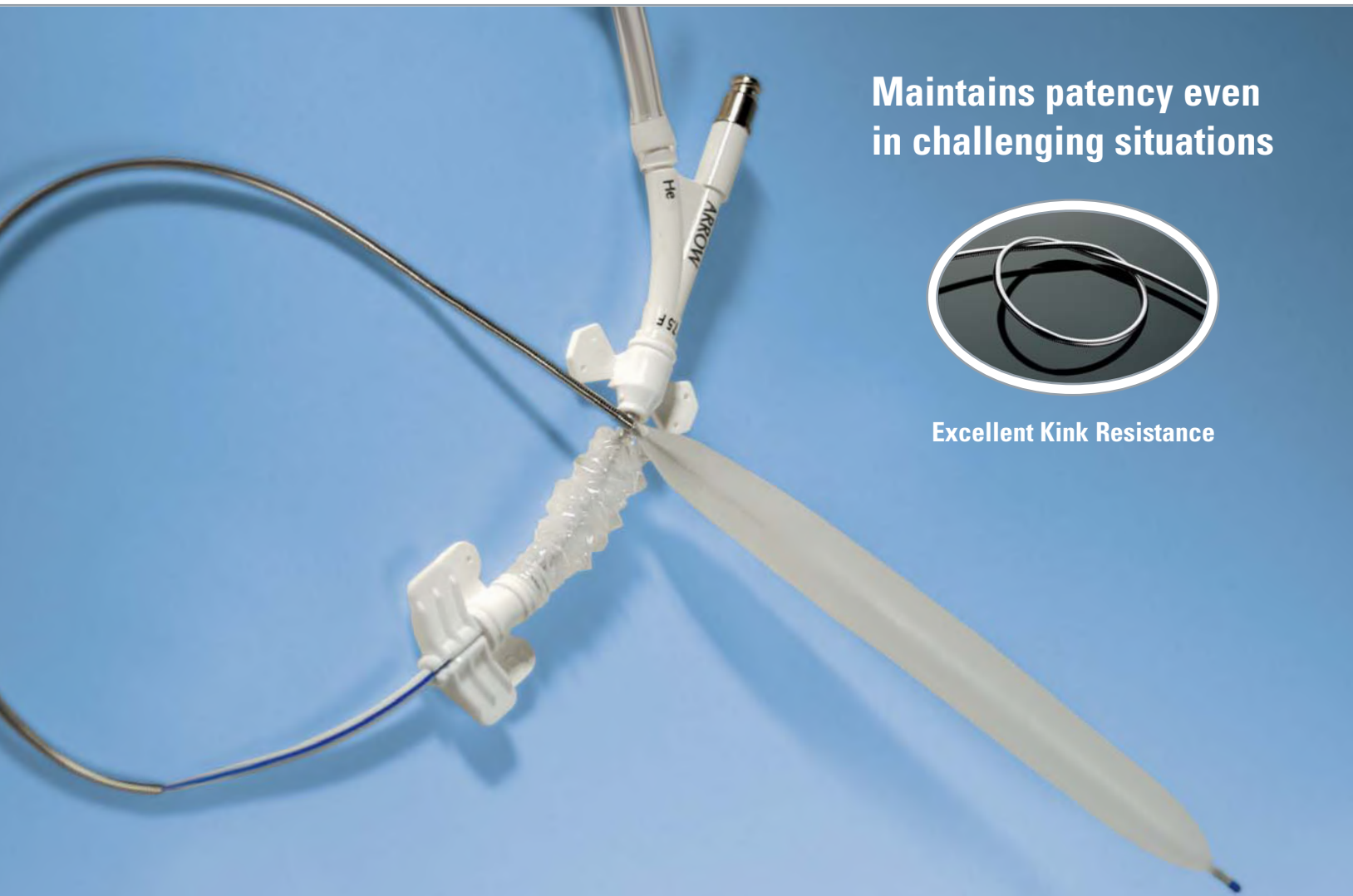

Maintains patency even
in challenging situations

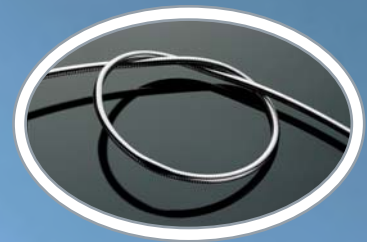

Excellent Kink Resistance


UltraFlex™ 7.5 Intra-Aortic Balloon Catheter

The Smallest Wire-Reinforced 30cc and 40cc IAB Catheter




Large Bore Connector

Maximizes gas flow.

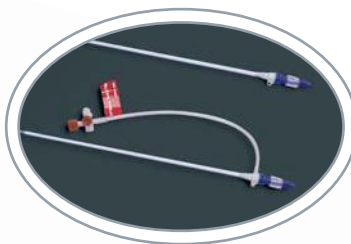


Optional Sheath with Sideport

Allows for arterial blood draws through sideport.

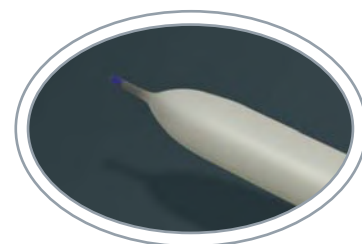


Large .027" Central Lumen

Minimizes clotting, reduces maintenance, and maximizes the quality of arterial pressure waveforms.

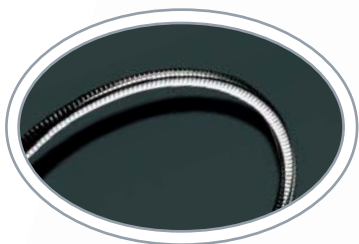


Kink Resistant Wire Reinforcement

Wire-reinforced catheter body will not collapse or kink, optimizing gas flow between pump console and balloon.

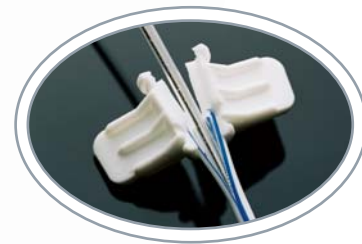


Optional Hemostasis Device

For post-insertion bleeding control.

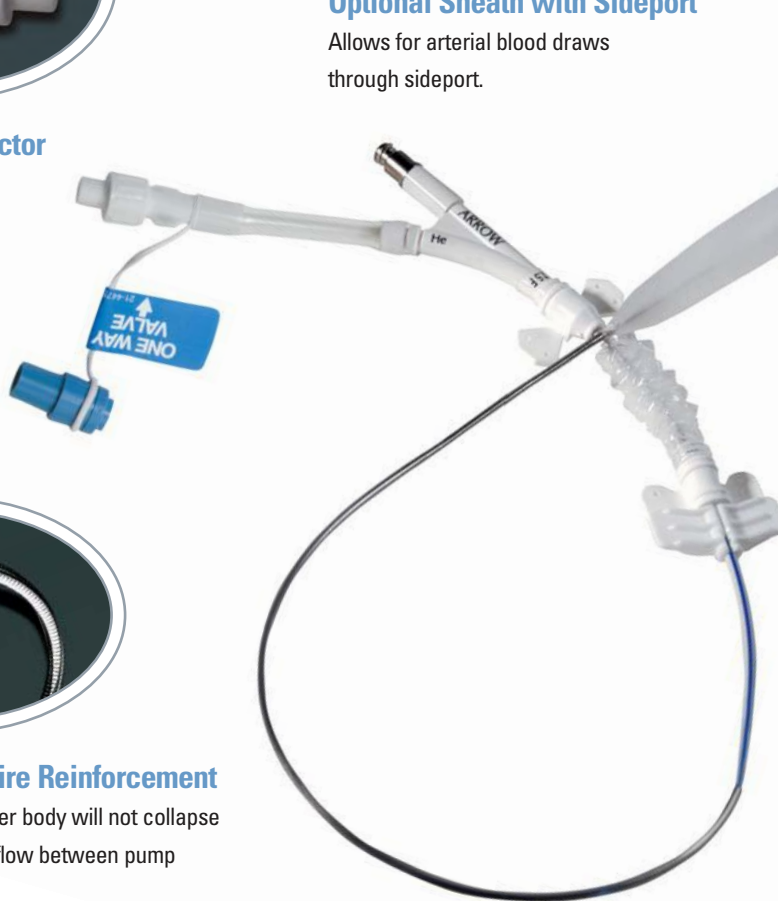


Wire-Reinforced and small 7.5 Fr. size delivers easy and reliable IAB performance

Kink Resistance Where You Need It Most

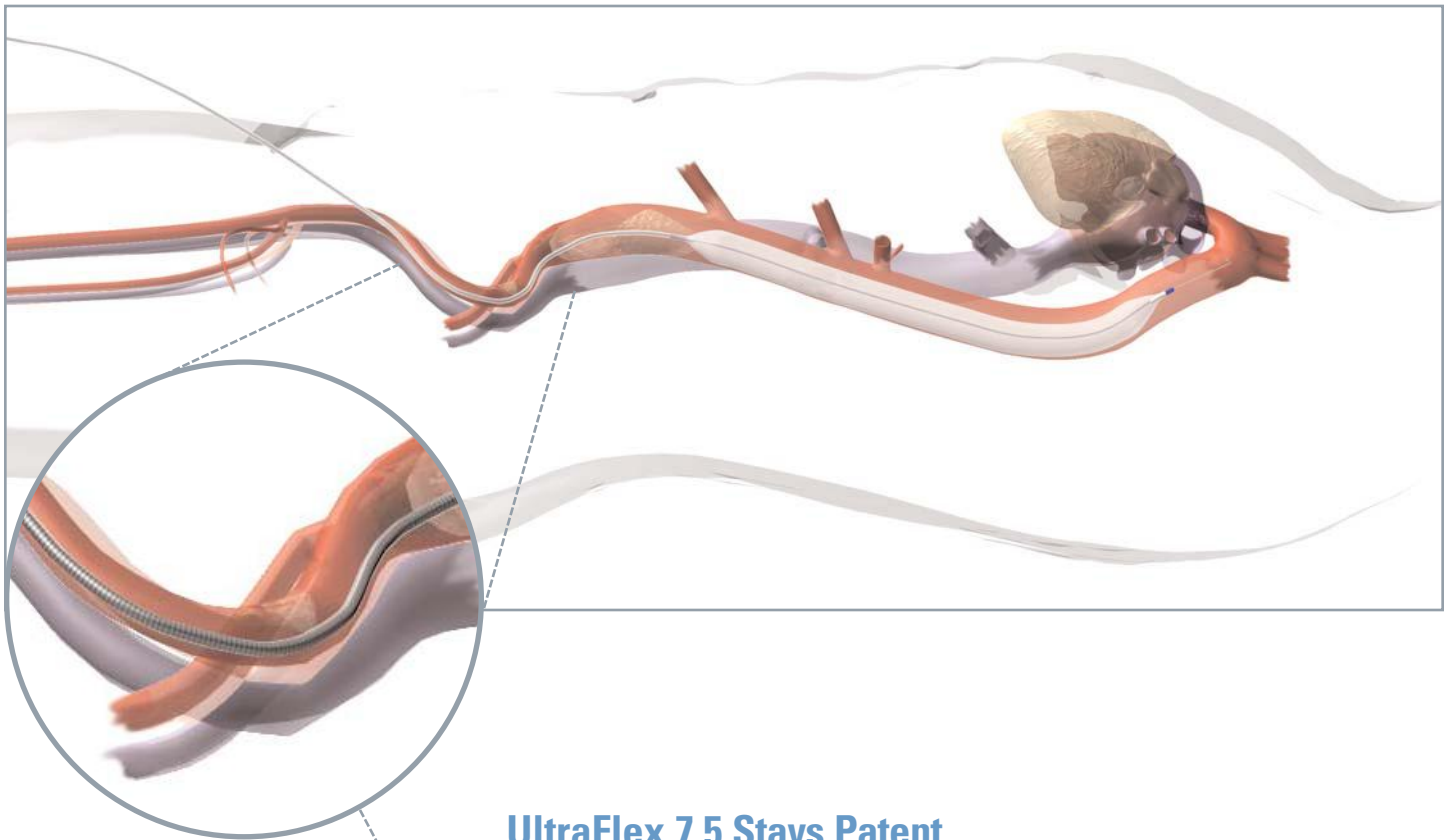

UltraFlex 7.5 Stays Patent

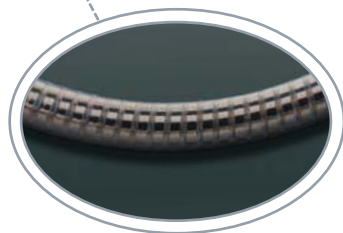

UltraFlex 7.5
Wire-reinforced UltraFlex 7.5
at a 20% bend.



Polyurethane catheter
Sample of competitor's
polyurethane catheter
at a 20% bend.

Product Specifications—UltraFlex 7.5 IAB Catheter

| Product Number | Catheter Size | Balloon Volume | Insertable Length Hemostasis Sheath | Catheter O.D. | Sheath Length | Central Lumen I.D. (inch) | Maximum Guidewire (inch) | Balloon Membrane Material | Catheter Material | Balloon Membrane Length | Inflated Diameter |
|----------------|---------------|----------------|-------------------------------------|---------------|---------------|---------------------------|--------------------------|---------------------------|-------------------|-------------------------|-------------------|
| IAB-06840-U | 7.5 Fr. | 40cc | 26.0" (66.0cm) | 7.5 Fr. | 6" (15cm) | .027 | .025 | Cardiothane™ II | SS/PU* | 10.2" (260mm) | 15mm |
| IAB-06830-U | 7.5 Fr. | 30cc | 23.8" (60.5cm) | 7.5 Fr. | 6" (15cm) | .027 | .025 | Cardiothane™ II | SS/PU* | 9.1" (230mm) | 13.9mm |
| IAB-06850-U | 9.0 Fr. | 50cc | 26.0" (66.0cm) | 9.0 Fr. | 6" (15cm) | .027 | .025 | Cardiothane™ II | SS/PU* | 10.3" (262mm) | 16mm |

ORDERING INFORMATION

Each balloon catheter is packaged with a pre-mounted hemostasis device, a 60cc syringe, one-way valve and driveline tubing for Arrow IABP consoles. Each IAB package contains an Insertion Kit (see components below) and separately packaged driveline tubing for TransAct® and Datascope® pump consoles.

| Product No. | Description |
|-------------|--|
| IAB-06840-U | UltraFlex 7.5 40cc Flexible Wire-Reinforced IAB Catheter. Pre-wrapped double-lumen 7.5 Fr. wire-reinforced IAB catheter with repositioning sleeve, pre-mounted "peel-away" hemostasis device for sheathless insertion. Insertion kit contains hemostasis sheath and hemostasis sheath with sideport for sheathed insertion. Driveline tubing with locking Quick-Connect®. |
| IAB-06830-U | UltraFlex 7.5 30cc Flexible Wire-Reinforced IAB Catheter. Pre-wrapped double-lumen 7.5 Fr. IAB wire-reinforced catheter with repositioning sleeve, pre-mounted "peel-away" hemostasis device for sheathless insertion. Insertion kit contains hemostasis sheath and hemostasis sheath with sideport for sheathed insertion. Driveline tubing with locking Quick-Connect®. |
| IAB-06850-U | UltraFlex 9 50cc Flexible Wire-Reinforced IAB Catheter. Pre-wrapped double-lumen 9.0 Fr. IAB wire-reinforced catheter with repositioning sleeve, pre-mounted "peel-away" hemostasis device for sheathless insertion. Insertion kit contains hemostasis sheath and hemostasis sheath with sideport for sheathed insertion. Driveline tubing with locking Quick-Connect®. |
| IAK-02691 | Driveline tubing with pre-attached Arrow pump connector for use with Arrow pump consoles and Arrow 30cc IAB catheters. |
| IAK-02692 | Driveline tubing with pre-attached Arrow pump connector for use with Arrow pump consoles and Arrow 40cc IAB catheters. |
| IAK-02693 | Driveline tubing with pre-attached Arrow pump connector for use with Arrow pump consoles and Arrow 50cc IAB catheters. |
| IAK-02262 | Driveline tubing for use with TransAct® pump consoles and Arrow IAB catheters. |
| IAK-02263 | Driveline tubing for use with Datascope balloon pump consoles and Arrow 30cc and 40cc IAB catheters. |
| IAK-02271 | Driveline tubing for use with Datascope balloon pump consoles and Arrow 50cc IAB catheters. |

REPLACEMENT INSERTION KITS

| | |
|------------------|--|
| IAK-06845 | Replacement Insertion Kit for use with UltraFlex 7.5, 30cc and 40cc (IAB-06830-U and IAB-06840-U) catheters. Each kit contains the following: <ul style="list-style-type: none"> • One: 18 Ga. x 2 1/2" Arterial Needle • One: 8.0 Fr. Sheath Dilator Assembly • One: 8.0 Fr. Sheath with Sideport and Dilator • Two: .025 x 175cm Teflon® Coated 3mm 'J' Extra Stiff Wire Guides • One: 8.0 Fr. Vessel Pre-Dilator • One: #11 Scalpel |
| IAK-06855 | Replacement Insertion Kit for use with UltraFlex 9, 50cc (IAB-06850-U) catheter. Each kit contains the following: <ul style="list-style-type: none"> • One: 18 Ga. x 2 1/2" Arterial Needle • One: 10 Fr. Sheath Dilator Assembly • One: 10 Fr. Sheath with Sideport and Dilator • Two: .025 x 175cm Teflon Coated 3mm 'J' Extra Stiff Wire Guides • One: 8.0 Fr. Vessel Pre-Dilator • One: #11 Scalpel |

*Wire-reinforced catheter body

Caution: Federal law restricts this device to sale by or on order of a physician. Contents of unopened, undamaged package are sterile. Disposable. Refer to package insert for current warnings, indications, contraindications, precautions, and instructions for use. Datascope® is a registered trademark of Datascope, Inc. Quick-Connect® is a registered trademark of Linvatec Corporation. Teflon® is a registered trademark of DuPont Corporation.



www.arrowintl.com

To locate an Arrow Distributor or for ordering information, contact Customer Service at:

United States

Tel: 800.523.8446

Fax: 800.343.2935

International

Tel: +1.610.655.8522

Fax: +1.610.655.8566

Distribution Worldwide:

Arrow offices are located in Belgium, Canada, Czech Republic, France, Germany, Greece, India, Italy, Japan, Mexico, Netherlands, Slovakia, South Africa, Spain, and the United States.

Corporate Headquarters address: Arrow International, Inc., 2400 Bernville Road, Reading, PA 19605 U.S.A. Tel: 610-378-0131

